


# ProSec® CIRCULAR STAPLER

ProSec® Circular Stapler is a single-use instrument that places two staggered circular rows of titanium staples and a circular knife divides the inner tissue. It is used in open and laparoscopic surgical procedures. Circular staplers are used for end-to-end anastomosis after bowel resection and in esophagogastric surgery.

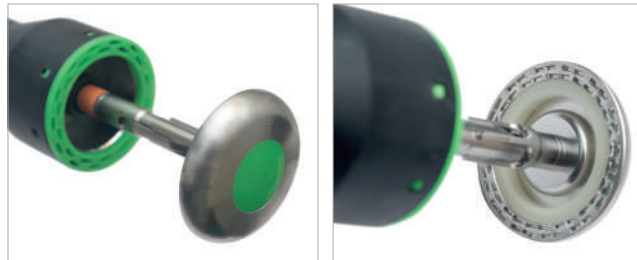

- ♦ Low profile anvil with atraumatic anvil edge for easy removal after firing
- ♦ Ring for audible click to confirm complete firing

## Features and Benefits

- ♦ Detachable anvil for ease of tying purse-string suture over shaft
- ♦ Extra staples per row for better anastomosis
- ♦ Extra long 51 cm shaft for versatile usage in the GI tract
- ♦ Airtight design for use in open and laparoscopic procedures
- ♦ Barlocked cartridge attachment for better connection and smooth placement
- ♦ Thin sloped anvil surface for easy removal

Ref. Code	Outer Diameter	Inner Diameter	Staple	Height	Formed Height	Colour
LPCS25X	25mm	17mm	22	4.8mm	2.0mm	●
LPCS29X	29mm	21mm	26	4.8mm	2.0mm	●
LPCS31X	31mm	23mm	28	5.0mm	2.2mm	●
LPCS33X	33mm	25mm	32	5.2mm	2.2mm	●