


ProSec® CURVED CUTTER

ProSec® Curved Cutter is a single-use instrument used in colo-rectal procedures that require deeper access within the pelvis. Its unique curved head places 2 double staggered rows of staples and a curved blade cuts in between. Prosec® Curved Cutter allows for placement of a 40mm cutline in the width of a 30mm space with low pelvic access.

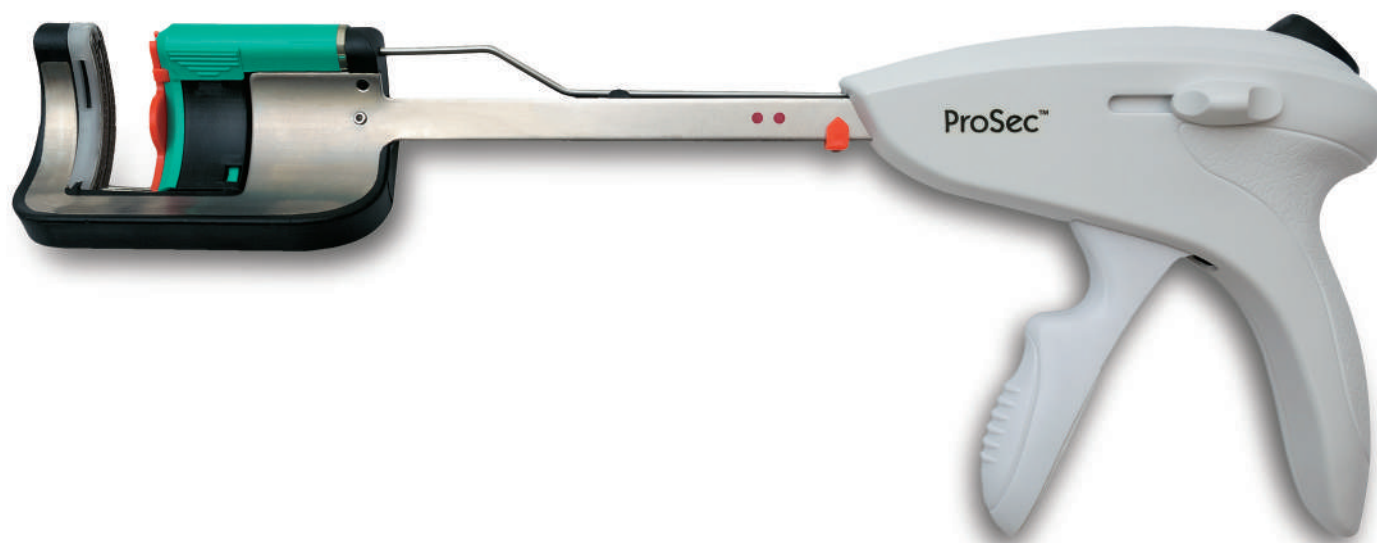

Features and Benefits

- ◆ Curved head enables deep pelvic access
- ◆ Easier placement for better tissue manipulation
- ◆ Ability to transect perpendicular to the rectum helps avoid a longer staple line
- ◆ Less tissue movement for more controlled transection
- ◆ Simultaneous stapling and cutting
- ◆ Optimized staple formation

Ref. Code	Description	Length	Staple Leg Length	Closed Staple Height
LPCC40G	Curved Cutter	40mm	4.8mm	2.0mm